


Kodak Medical X-Ray Film General Purpose Blue (MXB)

Get KODAK Quality in an Affordable “Blue” Film for General Radiography

Bring KODAK Film Technology to your x-ray department with a UV/blue-sensitive film that won't stretch your budget.

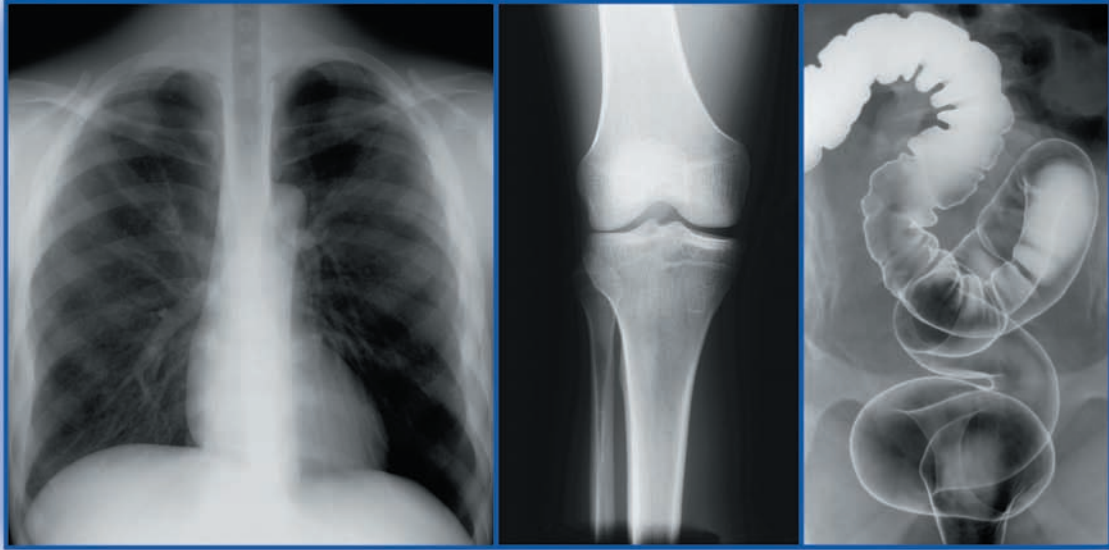

KODAK Medical X-Ray Film General Purpose Blue is an affordable, medium contrast film designed for general radiography.

Features include:

- KODAK Film quality and reliability
- Film speed among the highest of any blue-sensitive radiography film
- Enhanced blue image tone delivers optimal viewing characteristics and reduces eye fatigue
- Manual or automatic processing in standard cycle
- Excellent speed and contrast stability under variable processing conditions helps to ensure consistent results


Quality image-KODAK Technology

KODAK Medical X-Ray Film General Purpose Blue is part of a family of affordable Carestream Health offerings for radiography including films, screens, cassettes, film processing chemicals, and automatic film processors.



This UV/blue-sensitive film provides quality and reliability based on Kodak film technology and imaging expertise. Robust processing tolerance helps to deliver consistent results. Film speed, among the highest available with any blue-sensitive film, allows reduced exposure dose.

KODAK Medical X-Ray Film General Purpose Blue is designed for use with economical KODAK Medical X-Ray Screen/Blue 200. It is also designed to work well with a wide range of blue- and UV-emitting screens from other manufacturers.


Kodak Medical X-Ray Film General Purpose Blue (MXB)

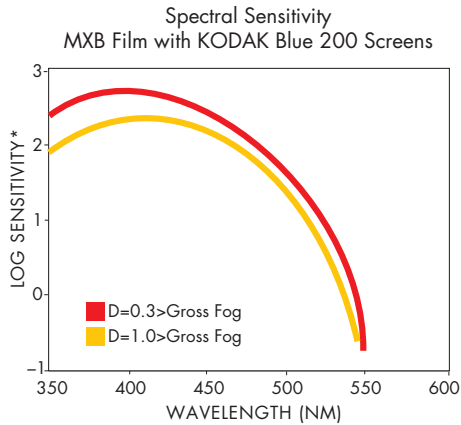

* Sensitivity=reciprocal of exposure (ergs/sq cm) required to produce specified density

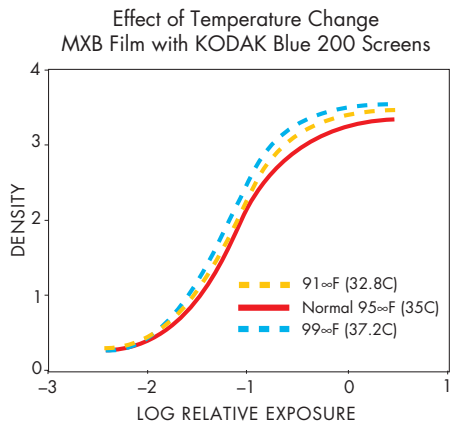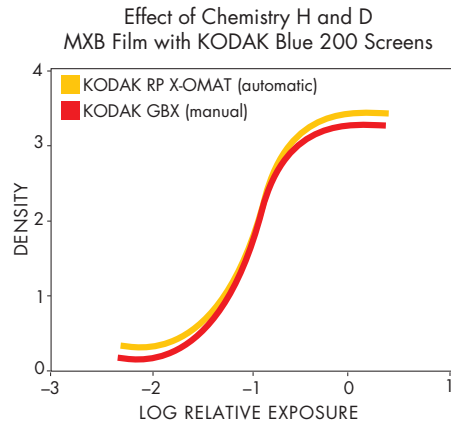

Film speed and contrast remain very stable across a wide range of developer temperatures

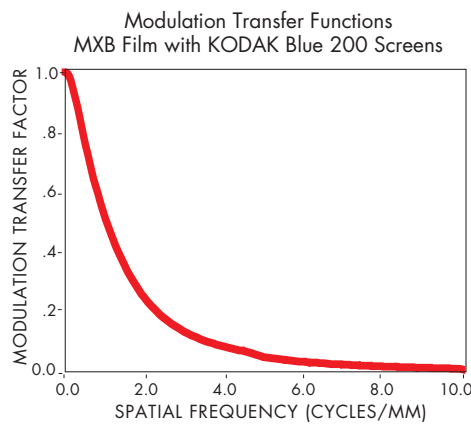

Specifications

Film Sizes	CAT No.
100 SH 35 x 43 cm	152 2200
100 SH 35 x 35 cm	172 8609
100 SH 30 x 40 cm	175 7699
100 SH 30 x 35 cm	850 5521
100 SH 24 x 30 cm	166 8144
100 SH 24 x 24 cm	890 3510
100 SH 18 x 24 cm	814 3059
100 SH 18 x 43 cm	175 7921
100 SH 11 x 14 in	181 8236
100 SH 10 x 12 in	820 5155
100 SH 8 x 10 in	822 5526
25 MXG-3 (tri-fold) SH 14 x 36 in	822 1746
25 MXG-3 (tri-fold) SH 14 x 51 in	869 0869

* See x-ray film distributor for current size offering

Film Processing

Manual processing or automatic processing using standard and rapid (60 second) cycles.

Film care and handling

Storage:

50 to 75° F (10 to 24° C),
30 to 50% relative humidity

Safelight:

KODAK GBX-2 or equivalent, standard frosted 15-watt bulb at least 4 feet (1.2 meters) from film-handling surface/area. The latest in darkroom lighting technology, the KODAK LED safelight (CAT. No. 852 6444) may be used. It provides twice as much visible light as conventional systems without compromising image quality.

Kodak

Licensed Product

Carestream Health, Inc.
150 Verona Street
Rochester, N.Y. 14608
U.S.A.

Carestream Health Canada Company
6 Monogram Place
Suite 200
Toronto, Ontario M9R 0A1
CANADA

T-Grain and X-Omatic are trademarks of Carestream Health, Inc.
The Kodak trademark and trade dress are used under license from Kodak.
Printed on recycled paper containing 10% post consumer waste fiber, using soybean-based inks.
M6-203 Printed in U.S.A. © Carestream Health, Inc. 2009 CAT No. 190 3442

For More Information

To learn more contact your Carestream Health representative or call 877-865-6325, Ext 538.

→ www.carestreamhealth.com

Carestream 
HEALTH